

Instructions For Use MagnetOs Putty synthetic bone void filler

DESCRIPTION

MagnetOs Putty is a synthetic 65-75% Tri-Calcium Phosphate (TCP - $\text{Ca}_3(\text{PO}_4)_2$) and 25-35% Hydroxyapatite (HA - $\text{Ca}_{10}(\text{PO}_4)_6(\text{OH})_2$) resorbable micro-structured bone void filler for the repair of bony defects. The product consists of granules, premixed with a synthetic polymeric binder that provides cohesion between the granules.

MagnetOs Putty is both osteoconductive and osteoinductive and has a porous trabecular structure that resembles the interconnected porosity of human cancellous bone. While the polymeric binder is rapidly resorbed after implantation, the granules of MagnetOs Putty induce and guide the three-dimensional regeneration of bone in the defect site into which it is implanted. When placed next to viable host bone, new bone will be deposited on the surface of the implant. The implant resorbs and is replaced by bone during the natural process of bone remodelling.

MagnetOs Putty is a ready for use product. Pressure applied by manipulation allows the shaping of MagnetOs Putty to conform to the defect contours.

MagnetOs Putty is gamma-sterilised, comes in several sizes in block form and is sterile packaged for single use only.

INTENDED USE

MagnetOs Putty is intended for use as bone void filler for voids and gaps that are not intrinsic to the stability of the bony structure. MagnetOs Putty is indicated for use in the treatment of surgically created osseous defects or osseous defects resulting from traumatic injury to the bone. MagnetOs Putty is intended to be packed into bony voids or gaps of the skeletal system (i.e. extremities, spine, cranial, mandible, maxilla and pelvis) and may be combined with autogenous bone.

MagnetOs Putty should not be used to treat large defects that in the surgeon's opinion would fail to heal spontaneously. In load bearing situations, MagnetOs Putty is to be used in conjunction with internal or external fixation devices.

Application in general bone surgery (entire skeletal system)

Replacing or supplementing autogenous / allogeneous spongiosa, e.g. for:

- Filling and bridging of skeletal bone defects including spine;
- Plastic reconstruction of damaged or resected bone areas;
- Filling of intervertebral implants.

Application in oral and maxillofacial surgery and dentistry

Filling or reconstruction of multiple walled (artificial or degenerative) bone defects, e.g.:

- Defects after the extirpation of bone cyst;
- Augmentation of an atrophied alveolar ridge;
- Sinus lift or sinus floor elevation;
- Filling of alveolar defects after tooth extraction for preservation of the alveolar ridge;
- Filling of extraction defects for creating an implant bed;
- Filling of two- or multiple-walled bone pockets as well as the bi- and trifurcations of teeth;
- Defects after operative removal of retained teeth or corrective osteotomies;
- Other multiple-walled bone defects of the alveolar processes and the facial skull.

CONTRAINDICATIONS

Use of MagnetOs Putty synthetic bone void filler is CONTRAINDICATED in the presence of one or more of the following clinical situations:

- To treat conditions in which general bone grafting is not advisable;
- In conditions where the surgical site may be subjected to excessive impact or stresses, including those beyond the load strength of fixation hardware (e.g. defect site stabilisation is not possible);
- In case of significant vascular impairment proximal to the graft site;
- In case of severe metabolic or systemic bone disorders that affect bone or wound healing;
- In case of acute and chronic infections in the operated area (soft tissue infections; inflamed, bacterial bone diseases; osteomyelitis);
- When intraoperative soft tissue coverage is not planned or possible;
- In direct contact with the articular space;
- In case of treatment with pharmaceuticals interfering with the calcium metabolism.

Despite the presence of some of the listed circumstances, the use of MagnetOs Putty may be the best solution for rectifying bone defects. The patient must be duly informed of the possible effects of complicating circumstances on the anticipated success of using MagnetOs Putty.

WARNINGS, CAUTIONS and PRECAUTIONS

Warning: MagnetOs Putty does not possess sufficient mechanical strength to support reduction of the defect site. Rigid fixation techniques are recommended as needed to assure stabilisation of the defect in all planes. MagnetOs Putty cannot be used to obtain purchase for screws. Screws must gain purchase in the host bone.

Warning: The granules contained in MagnetOs Putty must not be damaged or altered (e.g. by excessive compaction or crushing of the implant). Avoid overfilling of the defect as tension free wound closure is required.

Warning: As with any major surgical procedures, there are risks involved in bone grafting surgery such pain, swelling, superficial wound infection, deep wound infection, deep wound infection with osteomyelitis, loss of reduction, loss of bone graft, graft protrusion and/or dislodgment, and general complications associated with anaesthesia use and/or surgery.

Warning: Implantation of foreign materials can result in an inflammatory response or allergic reaction.

Caution: MagnetOs Putty is to be stored at ambient temperature (max 45 °C). Higher temperatures may affect the consistency and the ability of the device to retain its shape.

Caution: MagnetOs Putty's radiopacity is comparable to that of bone and diminishes as it is resorbed. This moderate radiopacity may mask underlying pathological conditions and must be considered when evaluating X-rays.

Caution: Inspect all packaging and components for damage before use. Do not use the device if it is damaged in any way.

Caution: Dosage is for SINGLE USE ONLY. Remaining material must be discarded. DO NOT re-use or re-sterilise. Re-use or re-sterilisation of the device may pose additional risks including, but not limited to, transmission of infectious agents.

Caution: Confirm expiration date before use. Do not use if expiration date has been exceeded.

DIRECTIONS FOR USE

1. The device can be shaped by finger manipulation to fit the defect contours.
2. The desired consistency and malleability can be achieved by pressure and warming in surgeon's hands overtime.
3. The product is ready for use: mixing with aqueous solutions is not recommended.
4. Defect site should be filled completely and secured to prevent migration of the implant.

MagnetOs Putty is intended for use by surgeons familiar with bone grafting and rigid fixation techniques. Familiarisation with the device and proper knowledge of bone grafting and rigid fixation techniques are extremely important.

Radiographic evaluation of the defect site is essential to accurately assess the extent of a traumatic defect and to aid in the selection and placement of the bone void filler and fixation devices. MagnetOs Putty must only be employed by or under the supervision of medical professionals with experience in the required surgical techniques and the use of biomaterials.

The exact operating procedures depend on the location, type and size of the defect. Close contact with vital bone is important for its function as a bone regeneration material and, therefore, a thorough freshening of the bone surface before applying the device is recommended (e.g. removal of bone fragments and necrotic tissue).

The defect must be completely filled with the device. Strong compacting or destruction of the granules structure (e.g. by crushing) contained in MagnetOs Putty must be avoided. Overfilling must be avoided to achieve a tension free closure. Fixation of the implant site must be sufficient to prevent collapse and deformity secondary to functional loading.

Anatomical reduction and rigid fixation in all planes must be obtained to ensure that the graft is not supporting load.

For large defects MagnetOs Putty may be admixed with spongiosa of comparable size. The selection of MagnetOs Putty size depends on the size of the defect to be filled.

Postoperative patient management should follow the same regimen as similar cases utilising autogenous bone grafting. Standard postoperative practices should be followed, particularly as applicable to defect repairs involving the use of fixation devices.

Particulars on application in oral and maxillofacial surgery and dentistry

For endosseous dental implants a time interval of 4 – 6 months should pass between defects filling with MagnetOs Putty and placing of the implant, in case of a sinus lift 6 months, if judged appropriate by the surgeon. In cases of larger defects surfaces in oral and maxillofacial surgery and dentistry, the user must decide on the use of membrane technique (GBR = Guided Bone Regeneration).

STERILISATION

MagnetOs Putty is provided sterile (Gamma irradiation). Do not re-sterilise.


HOW SUPPLIED

MagnetOs Putty is provided as a sterile, single use device. Do not use if package is opened or damaged.

Kuros Biosciences B.V.
Prof. Bronkhorstlaan 10, building 48
3723 MB Bilthoven
The Netherlands
www.kurosbiosciences.com
customerservice@kurosbio.com



Last revision of this text: 22Jan19

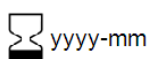

917-019 [04]